



ADEC MARINE LIMITED
4 MASONS AVENUE
CROYDON, SURREY, CRO 9XS
TEL 020 8686 9717

MOBMD  **T**

The logo icon for MOBMD is a stylized orange triangle pointing to the right, with a small gap between the 'D' and the 'T'.

Deployment instructions
Packing instructions

MOBMDT


Deployment instructions




1 Use available halyard or block and tackle to shackle both D rings



2 Mobmat being removed ready for deployment (see packing instructions)

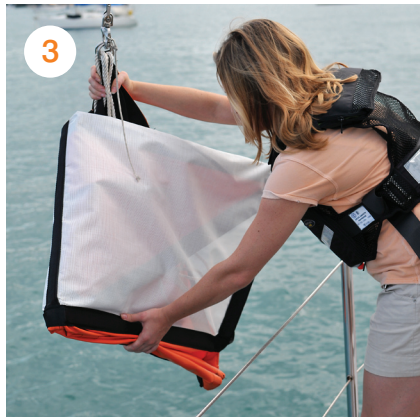

3 Deploy mobmat over the side

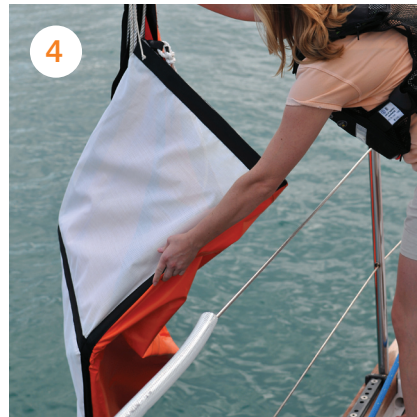


4 Solid orange PVC must be on the seaward side of the vessel

MOBMDT


Deployment instructions



5 Mobmat deploys effortlessly allowing the casualty to be recovered in a horizontal position

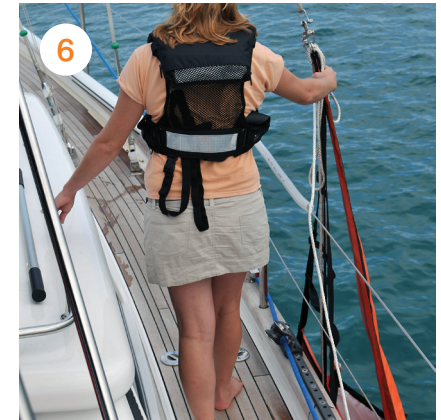

6 Mobmat can be repositioned if required

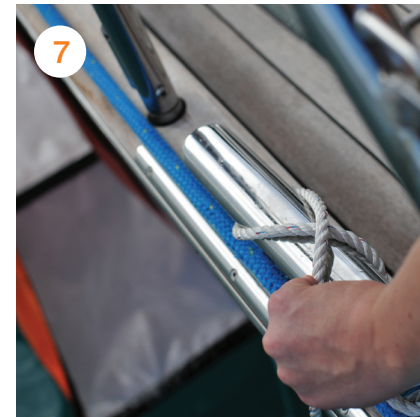

7 Attach provided lanyards to cleats or stanchions to steady Mobmat especially in high winds and seas



8 Mobmat with casualty can be recovered using a winch

MOBMDT

Deployment instructions




Alternatively Mobmat and casualty can be hoisted using a block and tackle arrangement



This picture shows Mobmat being fed underneath the top guardrail



Removing the bottom guardrail allows the casualty to be lifted onto deck more easily

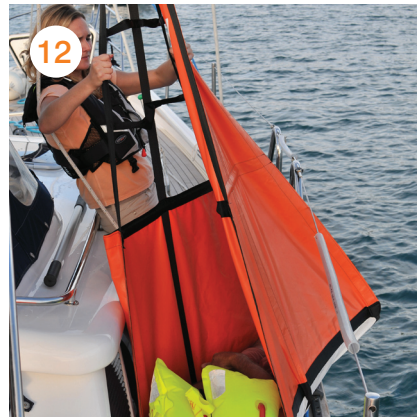

Mobmat recovers casualty comfortably especially if injured and reduces the risk of sudden shock and heart attack

MOBMDT


Packing instructions



Mobmat should be salt free and dry before packing



Note the lanyards already fixed to loops provided

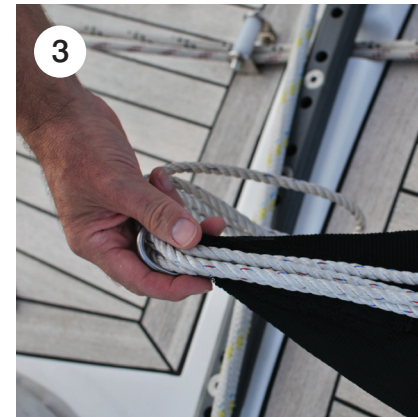

Lanyards loosely pulled through D rings ready for deployment



Take side battens and fold diagonally over one half of mesh base

MOBMDT

Packing instructions



5
Finish fold ensuring the battens are placed diagonally across one half of the mesh



6
Repeat with opposite side



7
Ensure all lanyards and webbing are flaked over the base



8
Grasp both D rings



9
Move D rings to a position ensuring they protrude for easy access



10
Fold base in half



11
Compress the two halves



12
Mobmat ready for deployment

MOBMDT

Packing instructions